
Aprilla

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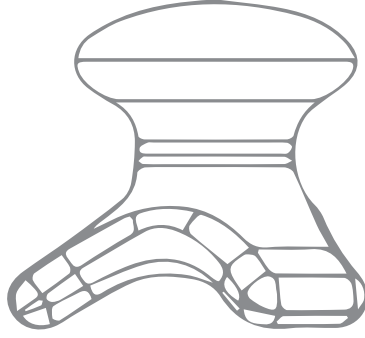TR EN

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Children Being Supervised not to Play with the Appliance

This appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliances by a person responsible for their safety. Children should be supervised to ensure that they do not play with the appliance.

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Hello,

First of all thank you for choosing our product. Please read this manual carefully before using your device. Please keep the manual as it may need in future. This guide is designed to help you to achieve maximum performance.

We wish you a pleasant use.


EN **Aprilla**

AFM 2803 User Manual

USING DEVICE

Battery Placement

- Open the top cover with care and insert 3 AAA thin pen batteries into the slot.
- Pay attention to the + and - directions when installing the batteries. Close the battery compartment

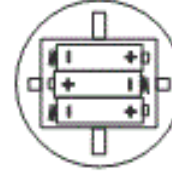Operating

- Bring the 3 legs of the massage device to the area where you want to massage and apply pressure. The vibration function will work automatically when you press the device.
- It will stop automatically when you stop applying pressure to the device.

FEATURES

- Three-point mini-massage tool ideal for neck pain, fibromyalgia pain, headache, back pain, legs, arms, face, shoulder pain and aching muscles.
- It will be better than any medication or ointment for your tense muscles.
- Very handy, small and lightweight.
- Portable, travel easy.
- One on / off switch design is easy to use, lightweight enough.
- Works with 3 AAA batteries (not included).

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
USING DEVICE

Battery Placement

- Open the top cover with care and insert 3 AAA thin pen batteries into the slot.
- Pay attention to the + and - directions when installing the batteries. Close the battery compartment.

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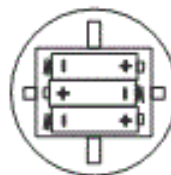FEATURES

Three-point mini massager perfect for neck pain, fibromyalgia pain, headache, back pain, legs, arms, face, shoulder pain and aches and pains. It will relax your muscles better than any medication or ointment. Very handy, portable, small and lightweight. One on / off switch, easy to use, it is sufficient to press slightly Works with 3 AAA batteries (not included).

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